

b.safe



HPLC All-In-One Box

A product brand of Bohlender

Strong NEWS


b.safe HPLC All-In-One Box

2

NOVELTIES

As HPLC already challenges lab experts, we help to simplify their work as much as possible. Our b.safe HPLC All-In-One Box includes all b.safe products that you need as basic equipment for one HPLC system.

A time-consuming search for suitable individual components is not necessary, and you can also postpone ordering the consumables (for the time being): Our set contains enough air valves and an exhaust air filter for a service life of 12 months.

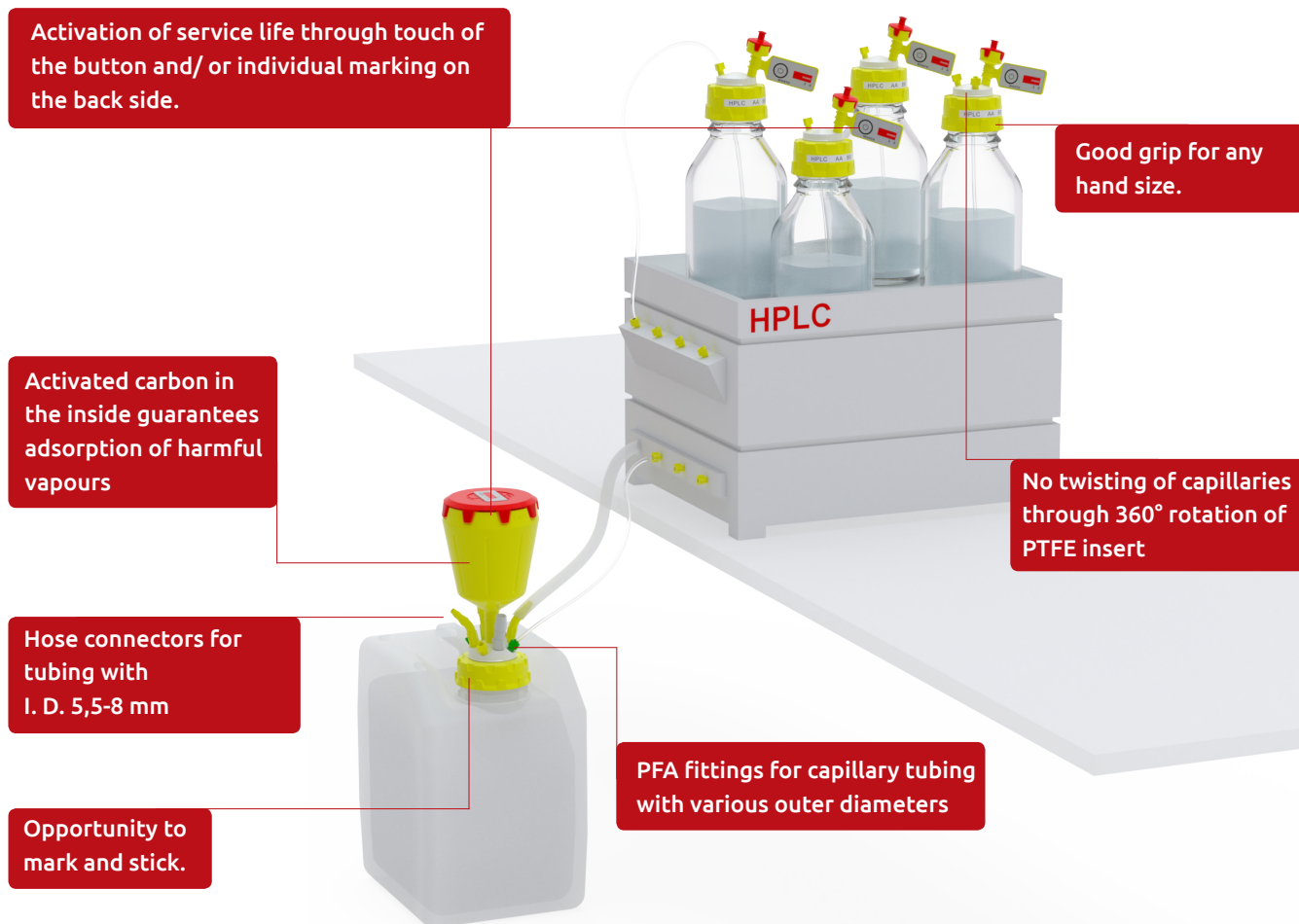

b.safe HPLC All-In-One Box

Material: **PTFE, PP** | Temperature resistance: **-20 °C to +200 °C** | chemical resistance: **+++ universal**

Complete set for equipping one HPLC system. Containing b.safe Caps GL45 for four laboratory flasks for solvent extraction as well as a b.safe Waste Cape with Exhaust Filter and a Canister for safe collection of solvent waste. The scope of delivery of the b.safe Caps includes PFA fittings for capillary O. D. 3,2 mm, Blind Fitting in PFA to close unused connections and Air Valves with service life indicator. For the collection of waste, the set includes a PE-canister S55 with 5-litre capacity, a suitable Waste Cap S55 and one Exhaust Filter L. The scope of delivery of the Waste Cap includes PFA Fittings for capillaries with OD 1,6/2,2/3,2 mm, Blind Fittings to close unused connections, tubing connectors for tubing I. D. 5,5-8 mm. Version and scope of delivery as per below chart. Consumables (Air Valves, Exhaust Filter) for a 12 months service life included.



Scope of delivery	Individual items for reordering:	Cat. No. All-In-One Box
A 3x b.safe Cap GL45 with 1x capillary connection, incl. 1x Fitting for Tubing O. D. 3,2 mm, 1 Air Valve with service life indicator.	M 145-01	R 9900-10
B 1x b.safe Cap GL45 with 2x capillary connection, incl. 2x Fitting for Tubing O. D. 3,2 mm, 1 Air Valve with service life indicator.	M 145-02	
C 4x b.safe Air Valve with service life indicator 6 months	M 505-01 (Packing unit 2 pcs.)	
D 1x b.safe Waste Cap S55 with 4x capillary connection, 4x hose connector, 1x connector for exhaust filter, incl. 4x Fitting for Tubing O. D. 1,6 mm (green) / 2,2 mm (purple) / 3,2 mm (yellow), 4x Blind Fitting for capillary connection, 4x Blind Fitting hose connection, 4x hose connector for tubing I. D. 5,5-8 mm	M 355-40	
E 1x b.safe Hose Connector for tubing I. D. 9-10,5 mm	M 503-16	
F 1x b.safe Exhaust Filter L , service life 12 months	M 506-03 (Packing unit 2 pcs.)	
G 1x b.safe Canister S55 , PE, capacity 5 litres	R 225-05	

Application:

Perfect set for the initial equipment of your HPLC system. Using b.safe Caps, eluents are directly extracted from laboratory flasks. The Air Valve prevents evaporation of solvents and only allows as much air to flow into the bottle as is necessary to equalize the pressure. Waste tubing are connected directly to the canister with the b.safe Waste Cap. The Exhaust Filter with activated carbon filling reliably restrains solvent vapours from the air escaping during filling.

b.safe

b.safe

A PRODUCT BRAND OF
BOHLENDER GmbH
Waltersberg 8
D 97947 Grünsfeld
Germany

Phone: +49 (0) 93 46 - 92 86-0
Mail: info@bohlender.de
www.bsafede.de

The current b.safe catalogue.

Huge choice on about 190 pages,
great service from the manufacturer.

Get your free catalogue!

www.bsafede.de/catalogue-request

b.safe

*Sophisticated
Safety Solutions
for HPLC.*

