

Non-Indicating Oxygen Traps

These oxygen/moisture traps are ideal for use in gas chromatography and GC/MS. These traps utilize a one-piece spun aluminium body and Swagelok™ end fittings. Each fitting incorporates a 40 micron stainless steel frit for particulate control.




	100cc	750cc
Oxygen Capacity	450mg	3,200mg
Oxygen Efficiency	<2 ppb	<2 ppb
Dimensions (L x D)	286 x 32mm	445 x 61mm
Pressure Rating	250 psi	250 psi

Description	Velocity Part Number		
	1/8" Fittings	1/4" Fittings	1/2" Fittings
Non-Indicating Oxygen Purifier, 100cc	TTO-100-2	TTO-100-4	-
High Capacity Oxygen Purifier, 750cc	TTO-750-2	TTO-750-4	TTO-750-8
Mounting Clips for 100cc Purifiers (2 required)	TTC-2		
Mounting Clips for 750cc Purifiers (2 required)	TTC-5		

Indicating Oxygen Traps

These glass-bodied traps utilize a nonhazardous material that reduces oxygen levels to less than 1 ppb. This indicating adsorbent changes colour from a light green to black as it adsorbs oxygen. Either of these traps is ideal as an indicator installed down stream of a Non-Indicating Oxygen trap. This allows high capacity removal of oxygen with the benefit of an indicator to alert the user of saturation of the primary trap. End fittings are nickel plated Swagelok™ with 40 micron stainless steel frits.

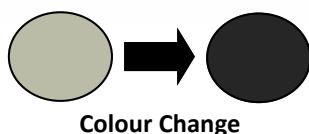

	50cc	150cc
Oxygen Capacity	50mg	150mg
Oxygen Efficiency	<1 ppb	<1 ppb
Dimensions (L x D)	242 x 29mm	242 x 38mm
Pressure Rating	100 psi	100 psi

Description	Velocity Part Number	
	1/8" Fittings	1/4" Fittings
Indicating Oxygen Purifiers, 50mg	TTIO-50-2	TTIO-50-4
Indicating Oxygen Purifiers, 150mg	TTIO-150-2	TTIO-150-4
Mounting Clips for 50cc Purifiers (2 required)	TTC-1	
Mounting Clips for 150cc Purifiers (2 required)	TTC-3	